

## **GLASS AGENCIES**

**5309, Anaj Mandi, Ambala Cantt. 133001**

**Haryana (INDIA)**

**Mobile : +91 98968-07858**

**E-mail : [glassagencies@gmail.com](mailto:glassagencies@gmail.com)**

**Website : [www.glassagencies.com](http://www.glassagencies.com)**



**Matrix**

Deals in Sales and  
Service of Ophthalmic  
Instruments

## **Auto Lensmeter Premio 30+**

**Matrix**


### **Salient Feature:**

- ▶ ISO standard 545mm green LED (108 points) with Hartmann sensor measures all types of lenses faster with more accuracy, consistency even without ABBE value compensation.
- ▶ Elegant design with improved quality & technology, ideal for versatile use.
- ▶ 7" Color TFT capacitive touch screen (820×480 resolutions)
- ▶ PD measurement, printer & UV measurement
- ▶ Auto detection of progressive power & Auto compensation of high-Index lenses
- ▶ Measurement with easy-refillable marking ink pad.


**M.: +91 9896807858**

## Premio 30+

The new advanced simultaneous measuring points with Hartmann sensor provide fast & reliable readings. The new advanced cluster of green LED measuring (545nm) beam.

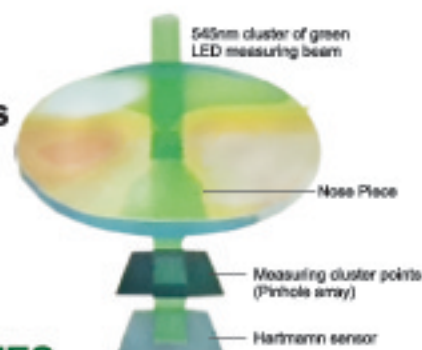

### PREMIUM TECHNOLOGY BEHIND PREMIO SERIES OF DIGITAL AUTO LENSMETER

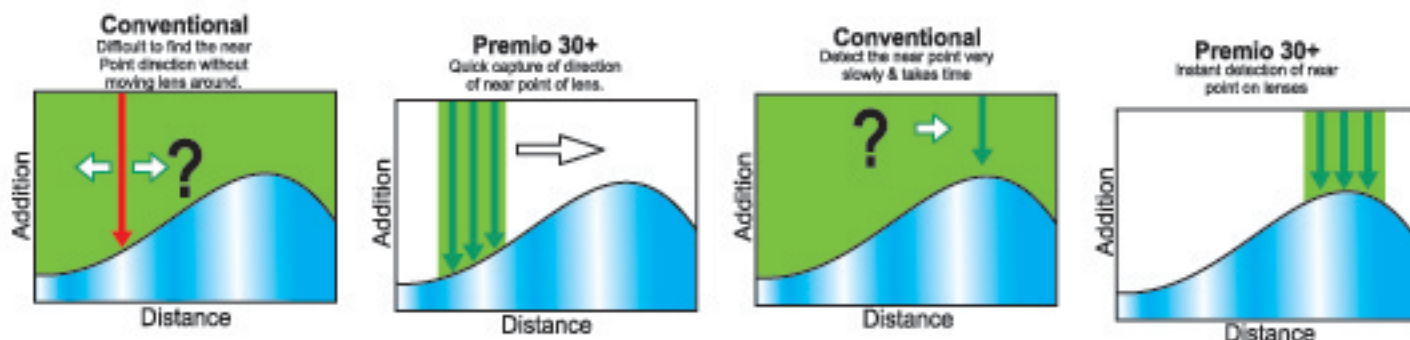

### Technical Specifications : Premio 30+

Measurement Range	Sph: $\pm 25.00D$ Cyl : $\pm 10D$ (Mode : +, -, mix) Step : 0.01D/0.06D/0.12D/0.25D Axis: 0~180° (step : 1°)
GLASS AGENCIES 5309, Anaj Mandi, Nr B. D. High School, Ambala Cantt. 133001 (Haryana) INDIA Mobile : +91 98968-07858	Addition : 0~+10D(step : 0.01/0.06/0.12/0.25D) Prism : 0 to 20 $\Delta$ Measuring Mode : Base, I/O,U/D (Cartesian co-ordinate) P-B : Polar co-ordinates (Step : 0.01 $\Delta$ )
Measurable Lens Dia	$\varnothing 20mm \sim \varnothing 120mm$ (Spec. Lens) Contact Lens : more than 5mm Dia
Measurable Lens	Frame Lens; Single Lens; Bifocal; Multifocal, Progressive Lens, ; Hard & Soft Contact Lens (with proper holder)
PD Measurement	42-82mm by easy scale cursor movement. P.H Range : 10mm~ 45mm
Display	7" color TFT touch screen
Printer	Thermal printer : 58mm width
ABBE compensation	Not required / Automatic
Interface Connector	RS 232/USB & ETHNET
Power Rating	100~240 VAC / 50~60Hz. Rated power : 35 VA
Dimensions/ Net wt.	H(44) W(20) D(25) cm/ 5kg
Packing/ Gross wt.	(L)55 $\times$ 38(W) $\times$ 33(H) cm/ 6.5kg
Environmental Condition	Temp : 0° C~40°C / Humidity : $\leq 85\%$