


# EROSE LIQUID LIMIT DEVICE

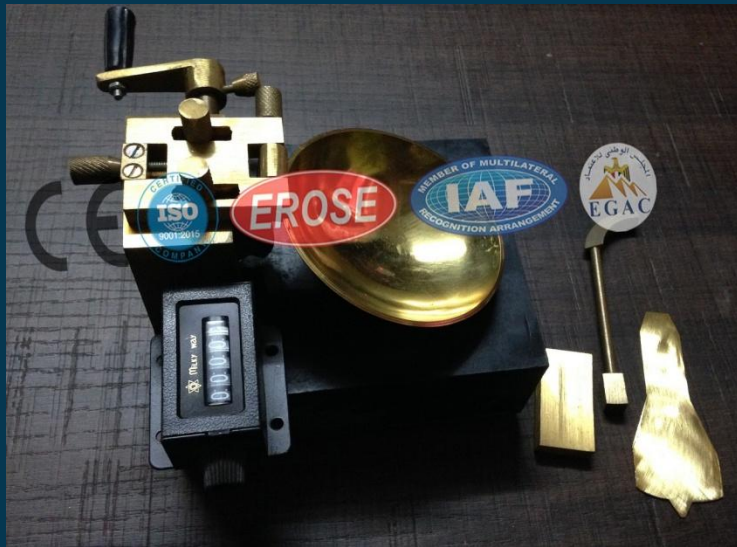

Glass Agencies are  
Manufacturer of  
**LIQUID LIMIT  
DEVICE** Under  
EROSE Brand

## SPECIFICATION

- Power Source : Hydraulic
- Material : SS, Brass
- Finishing : Polished
- Manufacturer Name : Glass Agencies
- Brand : EROSE
- Base : Moulded Hard Rubber
- Cup : Spun Brass, Pin & Slot Mounting
- Counter : Records number of cup Droops
- Color : Black and golden
- Use for : Soil Testing
- Feature : Easy to operate
- Weight : 2 kg

❖ **Liquid Limit Device:** fitted with counter complete with Casagrande grooving tools and a Gauge Block. The liquid limit is defined as the moisture content at which the soil passes from plastic state to liquid state; it is very helpful in the classification of the potential properties of the soil material. The liquid limit device consists of a hard rubber base carrying a sliding carriage assembly to which a brass cup is hinged. The cup is raised and allowed to fall through a height of 1cm on to the hard rubber base, by the help of a lead screw provided at the back of the sliding carriage.

## GLASS AGENCIES

5309, Anaj Mandi, Ambala Cantt. 133001 (Haryana) INDIA

Office : +91 89304-07858

Mobile : +91 98968-07858

E-mail : [glassagencies@gmail.com](mailto:glassagencies@gmail.com)

Website : [www.glassagencies.com](http://www.glassagencies.com)

